


EASYCLEAN




DRYING ASSEMBLY FOR SYRINGE SHIELDS WITH INTEGRATED DRIP TRAY

ALL ENERGIES

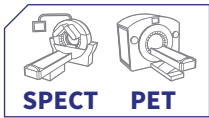

Easyclean is a new must-have accessory in the Lemer Pax nuclear medicine range. Focusing on **constant improvement of the daily applications** of our users, the support can be **directly immersed** in the soaking, disinfection and decontamination tanks of a minimum volume of 5 L. It is made up of 11 or 18 PVC rods depending on the model, each one associated with a shock absorbing washer, in order to avoid any risk of scratching and breaking of lead glass when inserting the syringe shield. It is possible to position respectively up to **18 syringe shields** simultaneously and thus avoid the usual risks of syringe shield breakage due to impact during soaking.

EFFECTIVE DIMENSIONS (mm)

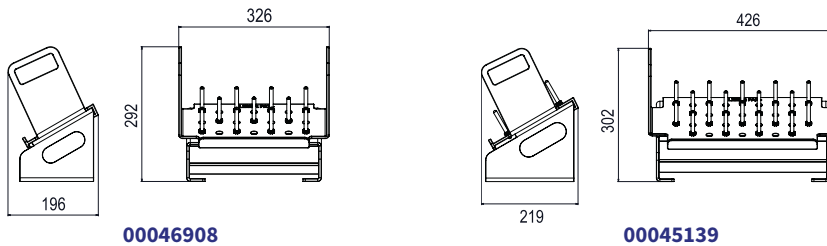

CHARACTERISTICS

General	Small model	Large model
Total dimensions of the assembly:	L 326 x D 196 x H 292 mm	L 426 x D 219 x H 302 mm
Weight of the assembly:	1,6 kg	2,3 kg
Drying rack		
Overall dimensions:	L 426 x P 220 x H 180 mm	L 426 x P 216 x H 187 mm
Materials:	White PVC	
Drying rods:	11 screwed in rods with rubber buffers dia. 8 mm	18 screwed in rods with rubber buffers dia. 8 mm
Weight:	1,17 kg	1,15 kg
Dryer support		
Overall dimensions:	L 292 x D 160 x H 158 mm	L 392 x D 200 x H 182 mm
Angle:	30°	
Materials:	White PVC	
Weight:	1,1 kg	1,14 kg
Ref. of the assembly:	00046908	00045139
Ref. drying rack:	00046906	00046901
Ref dryer support:	00047128	00045083
Ref. rods and rubber buffers:	00049081	

Drying syringe shields is made easy with a **benchtop drying rack**. The 5° angle for the dryer and 30° angle for its support allow for quick draining of the rinsing liquid, which is directly stored in the 1-litre retention tank. Each **lightweight** element features **two ergonomic handles**, allowing the syringe shield drying unit to be moved very easily. The choice of white PVC materials **simplifies daily maintenance**, and ensures an **extremely long service life**, without chemical interaction with hospital decontamination products or oxidization problems.

ASSOCIATED PRODUCTS

- Easyview syringe shields
- Mediclic syringe shields
- Easysqueeze syringe shields